



Finecare™ Full-automatic FIA Meter



Finecare™ Multi-channel FIA Meter

Powerful Features:

- High throughput:** 180 specimens/hour
- Full-automatic:** place the samples into the sample acceptor, click "test" and get results
- Mass memory:** 20G with result searching function
- Superior community:** LIS link, PC/USB interface, external printer, scanner (optional)

High throughput, Full-automatic, Continuous, Rapid, Accurate quantitative. Click and test, only 2 steps.

Powerful Features:

- Multi-channel:** 2-6 channels as required
- Multi tests:** 2-6 different test items at one time
- User-friendly interface:** 7 inch touch screen
- Superior communication:** LIS link, PC/USB interface, external printer, scanner (optional)

Multi-channel, multi-tests, rapid, full-scale quantitative test. Provide you with integrated diagnosis basis

High-standard quality system



Guangzhou Wondfo Biotech Co., Ltd.
No. 8 Lizhishan Road Scientific City, Luogang District, Guangzhou China P.R. 10663
Tel: +86-20-32296082 Fax: +86-20-32296062
E-mail: sales@wondfo.com.cn Websites: www.wondfo.com.cn

WONDFO BIOTECH
WeAreWorkingForYourHealth



Finecare™ FIA Meter

Fluorescence Immunoassay Rapid Quantitative Test system



Fineware™ FIA Meter

Fineware™ FIA Meter is a fluorescence immunochromatographic analysing system which can help diagnose conditions such as infection, diabetes, cardiovascular diseases, renal injury and cancers, etc.

The Fineware platform provides reliable and quantitative results of various kinds of analytes in human blood or urine within several minutes.

Application: Laboratory, ER, Cardiology, ICU, Respiratory, Pediatrics, etc.

Principle	Fluorescent immunochromatography
Physical Parameters	Dimensions: 250(L) x 258(W) x 104(H)mm Weight: 2.5kg Power adapter: 100-240AC, 50-60Hz
Optical description	Excitation light source: LED, 470nm Detector: 525nm Photodiode
I/O	LCD touch screen: 3.5" Printer: built-in thermal printer Information input: touch screen, barcode reader Record storage: 10,000 records Interface: RS232, USB, Ethernet (Direct connection to LIS / HIS)



Fineware™ Test Markers



Product Class	Test Item	Diseases	Specimen Type	Reaction Time	Linear Range	Clinical Reference	Qualification
Infection	PCT	Sepsis	WB/ Serum/ Plasma	15 min	0.1-100 ng/mL	0-0.5 ng/mL	CE CFDA
	CRP (hsCRP+CRP)	Inflammation	WB/ Serum/ Plasma	3 min	0.5-200 mg/L	CRP:0-10 mg/L hsCRP:0-1 mg/L	CE CFDA
Coagulation	D-Dimer	Thrombosis	WB/ Plasma	3 min	0.1-10 mg/L	0-0.5 mg/L	CE CFDA
Cardiac	NT-proBNP	Heart Failure	WB/ Serum/ Plasma	15 min	18-35000 ng/L	0-300 ng/L	CE CFDA
	cTnI	Myocardial Infarction	WB/ Serum/ Plasma	15 min	0.1-50 ng/mL	0-0.3 ng/mL	CE CFDA
	Myo	Myocardial Infarction	WB/ Serum/ Plasma	15 min	2.0-400 ng/mL	0-58 ng/mL	CE CFDA
	CK-MB	Myocardial Infarction	WB/ Serum/ Plasma	15 min	0.3-100 ng/mL	0-5 ng/mL	CE CFDA
	H-FABP	Myocardial Infarction	WB/ Serum/ Plasma	15 min	2-120 ng/mL	0-7 ng/mL	CE CFDA
	3 in 1 (NT-proBNP / cTnI /CK-MB)	Myocardial Infarction /Heart Failure	WB/ Serum/ Plasma	15 min	Same with single items	Same with single items	CFDA
Tumor	3 in 1 (cTnI/Myo/CK-MB)	Myocardial Infarction	WB/ Serum/ Plasma	15 min	Same with single items	Same with single items	CE CFDA
	AFP	Liver Cancer	Serum	15 min	5-400 ng/mL	0-20 ng/mL	CE CFDA
Diabetes & Renal Injury	PSA	Prostate Cancer	Serum	15 min	2-100 ng/mL	0-4 ng/mL	CFDA
	HbA1c	Diabetes Mellitus	WB	5 min	4.0-14.5%	0-6.5%	CE CFDA NGSP
	MAU	Renal Injury	Urine	10 min	5-300 ng/mL	0-20 ng/mL	CE CFDA
	Cys-C	Renal Injury	WB/ Serum/ Plasma	5 min	0.2-10 mg/L	0.5-1.1 mg/L	Applying

Choose Fineware™ For



Extra quick operating platform for large number of specimens



Powerful information management. (Direct connection to LIS)



Convenient and fast. Get results in 3-15 minutes

Unique quick mode

Designed sampler

Information management


Reliable accuracy

Instant result

Room temperature storage



Special designed sampler for small volume of specimen

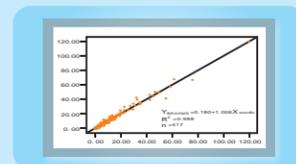

Good correlation with wet chemistry method results



Solid synthetic craft, can be stored at room temperature (4-30°C)

Principle and operation

